

Diaphragm seal MD 56 NEUMO BioControl®



Benefits

- Defined installation and press-fitting with metal stop
- Medium cannot get behind the seal
- SIP/CIP enabled

Application

Flange type diaphragm seal MD 56 for sterile process separation of pressure measuring and control units in pipes or tanks. For mounting to Bourdon tube pressure gauges, pressure transducers and pressure switches. For installation in sterile processes without dead spaces. Typical application areas:

Sterile applications in process engineering

Pharmaceutical industry

Biotechnology

Food and beverages industry

Beverage machines

Dairies

Breweries

Versions

NEUMO BioControl®

	Process connection	Instrument con- nection	Minimum range* (in bar) at nominal size				
			63	100	160	DMU	Part no.
	D25	Welded connection	4 bar	4 bar	0.6 bar	0.6 bar	31410W
	D50	Welded connection	0.6 bar	0.6 bar	0.6 bar	0.6 bar	31411W
	D65	Welded connection	0.6 bar	0.6 bar	0.6 bar	0.6 bar	31412W
	D80	Welded connection	0.6 bar	0.6 bar	0.6 bar	0.6 bar	31413W

Blue part no. = in-stock items





Description

All materials used are FDA-listed. Defined installation as per EHEDG hygienic design recommendations. The chemical seal and the measuring instrument are welded together, resulting in an inseparable single, shock- and vibration-resistant unit without external edges that might collect dirt. The chemical seals are available with nominal diameters 25, 50, 65, 80; due to their compact design, they are suitable for a wide range of applications.

Technical specifications

Process connection

Stainless steel 316 L

NEUMO BioControl®

D25, D50, D65, D80

Diaphragm

Stainless steel 316 L

Material (wetted parts)

Stainless steel 316 L, EPDM

USP CLASS VI classified, FDA-listed

Surface roughness RA

≤ 0.8 µm

Seal

EPDM

Options

- Other materials
- Electropolished
- Cooling element
- Capillary
- Other filling liquids
- Accessories such as housing, block flange and seals

Instrument connection

Welded connection

Pressure transmission liquid

Paraffin oil (FM 09), FDA-listed

Pressure range


See overview of versions

Can be mounted to


- Bourdon tube pressure gauges NG 50, 63, 80, 100, 160
- Pressure transducers (DMU)
- Pressure switches (DS)

Technical drawings

Pressure gauge NG 100 with MD 56

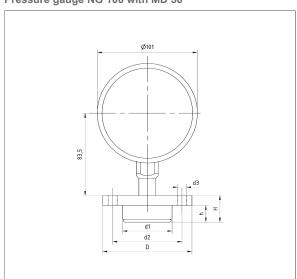

Pressure transducer DMU 02 Vario BC

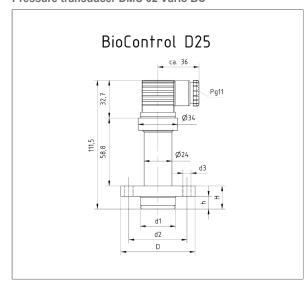





An- schluss	D	d1	d2	d3	Н	h
BioCon- troll® D25	64	30,4	50	4 x ∅ 7	20	11
BioCon- troll® D50	90	49,9	70	4 x ∅ 9	27	17
BioCon- troll® D65	120	67,9	95	4 x ∅ 11	27	17
BioCon- troll® D80	140	87,4	115	4 x ∅ 11	37	25

